




# DISPOSABLE SPREADERS

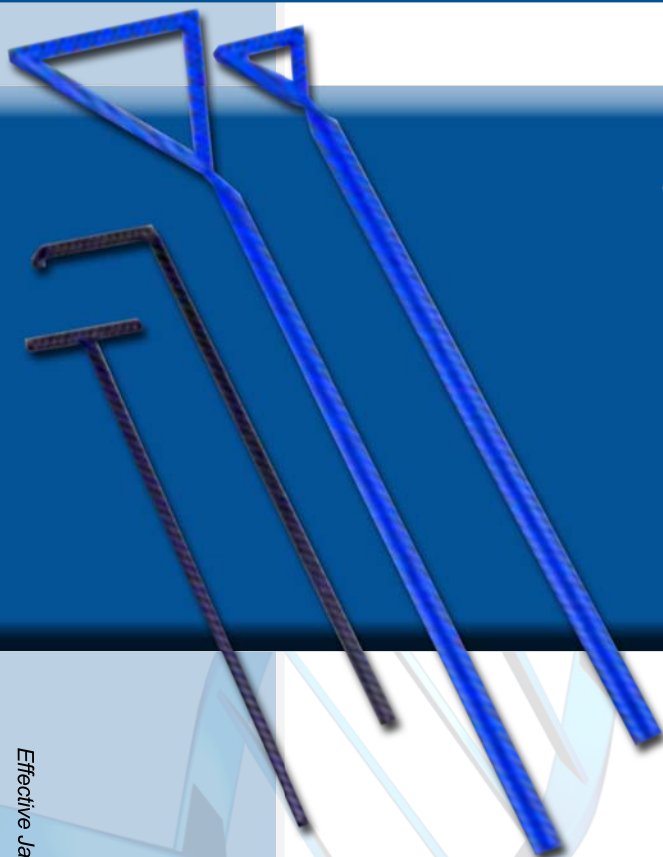

- **Spread cultures & liquids evenly**
- **Smooth edges - minimize "gouging"**
- **3 Configurations: "L", "T", or Triangle**
- **Ideal for bacterial CFU counts**
- **Presterilized - no flaming required**

The MSP line of bacteriological spreaders has been expanded to include four different shapes and sizes. In addition to the original L-shape, a similar T-shape version and two triangle shapes (large and small) are now available. All versions are supplied in sterile packaging, but are also autoclavable and have smooth, rounded edges to prevent unwanted cutting, scraping or gouging of agar during inoculation.

**L-shape** The original spreaders are ideal for spreading bacteriological cultures evenly. The shape allows the user to rotate and apply easily throughout an entire round plate or dish. They are made of ABS and have a slight upward curve at the edge to prevent dragging or scraping.

**T-shape** This shape ensures that consistently even pressure can be applied across the entire span of the 34mm spreading surface. Designed for spreading and dispersing liquids across the surface of agar culture plates, they are made of ABS and commonly used for CFU plate counts on water, milk, urine, etc.

**Triangle Shape** This new design is also used for CFU count applications, but has a grooved handle that improves grip and orientation. Made of polypropylene, they are available with a 30mm span or an extra large 60mm span.

## MedSupply Partners

3715 Atlanta Industrial Pkwy, Ste B  
Atlanta, Georgia 30331  
Phone: 1-877-633-7779  
[www.medsupplypartners.com](http://www.medsupplypartners.com)

### Ordering Information

Cat. No.	Description	Dimensions	Sterile Unit	Units/Pack	Qty/Case	Price
62-1014-1	Spreader, L-shape	38 x 156mm	10 per peel pack	50 peels	500	\$77.00
62-1014-2	Spreader, L-shape	38 x 156mm	1 per peel	500 peels	500	\$89.00
62-1014-3	Spreader, T-shape	34 x 140mm	10 per peel pack	50 peels	500	\$85.00
62-1014-4	Spreader, T-shape	34 x 140mm	1 per peel	500 peels	500	\$95.00
62-1014-5	Spreader, Triangle	60 x 235mm	25 per peel pack	4 peels	100	\$59.00
62-1014-6	Spreader, Triangle	60 x 235mm	1 per peel	500 peels	500	\$320.00
62-1014-7	Spreader, Triangle	30 x 207mm	25 per peel pack	4 peels	100	\$59.00